

GENIE® MATERIAL

GENIE® VPS MATERIAL 2-PACKS

2x 50mL Cartridges; 6 Mixing tips
(XLB & LB varieties include 6 intraoral tips)

Viscosity	Rapid Set (2:10)	Standard Set (4:30)
Extra Light Body	REF 77600	REF 77605
Light Body	REF 77610	REF 77615
Regular Body	REF 77620	REF 77625
Heavy Body	REF 77630	REF 77635

GENIE® VPS MATERIAL BULK KITS

60x 50mL Cartridges

Viscosity	Rapid Set (2:10)	Standard Set (4:30)
Light Body	REF 77720	REF 77725
Regular Body	REF 77700	REF 77705
Heavy Body	REF 77710	REF 77715

GENIE® MAGIC MIX

1x 380mL Cartridge; 10 mixing tips

Viscosity	Rapid Set (2:10)	Standard Set (4:30)
Regular Body	REF 78800	REF 78805
Heavy Body	REF 78810	REF 78815

VERSA-TEMP® 2 TEMPORARY MATERIAL

VERSA-TEMP® 2 TEMPORARY CROWN & BRIDGE RESIN

1 Cartridge Pack 1x 50mL Cartridge; 15 Mixing Tips	Shade	3 Cartridge Pack 3x 50mL Cartridges
REF 72041	BL	
REF 72051	A1	REF 72066
REF 72056	A2	REF 72071
REF 72061	A3	REF 72076

SILGIMIX® MATERIAL

SILGIMIX® ALGINATE REPLACEMENT MATERIAL

380mL Single Pack 1x 380mL Cartridge; 10 mixing tips	REF 79010
380mL Bulk Pack 4x 380mL Cartridges	REF 79015
50mL 8-Pack 8x 50mL Cartridges	REF 79020
50mL Bulk Pack 60x 50mL Cartridges	REF 79025

GENIE® MAGIC MIX BULK KITS

4x 380mL Cartridges

Regular Body	REF 78820	REF 78825
Heavy Body	REF 78830	REF 78835

GENIE® BITE MATERIAL

2x 50mL Cartridges; 6 blue mixing tips

	Super Fast Set (0:60)	Fast Set (0:90)
Bite Material	REF 77645	REF 77640

GENIE® PUTTY

2x 300mL jars; 2 scoops

	Rapid Set (2:30)	Standard Set (4:30)
PUTTY	REF 77650	REF 77655

IMPRESSION MATERIAL ACCESSORIES

BLUE Mixing Tips (50 ct)	REF 0072100FG
RED Mixing Tips (50 ct)	REF 0072101FG
RED Intraoral Mixing Tips (100 ct)	REF 0072008FG
Dynamic Mixing Tips (50 ct)	REF 0072012FG
50mL Dispenser Gun (S-50)	REF 70051
Tray Adhesive	REF 77100

SENSITEMP® NEZO TEMPORARY CEMENT

NON-EUGENOL ZINC OXIDE TEMPORARY CEMENT

5mL Syringe Pack 1x 5mL syringe; 10 mixing tips;	REF 70050
25mL Cartridge Pack 1x 25mL Cartridge; 15 mixing tips;	REF 70055


SENSITEMP® RESIN TEMPORARY CEMENT

NON-EUGENOL RESIN TEMPORARY CEMENT

5mL Syringe Pack 1x 5mL syringe; 10 mixing tips;	REF 70010
25mL Cartridge Pack 1x 25mL Cartridge; 15 mixing tips;	REF 70020

TEMPORARY CEMENT ACCESSORIES

20x Blunt Mixing Tips (Brown)	REF 70030
20x Needle Mixing Tips (Brown)	REF 70031
50x Mixing Tips (Blue) for 25mL Cartridge	REF 70035
25mL Dispenser Gun (S-25)	REF 70025

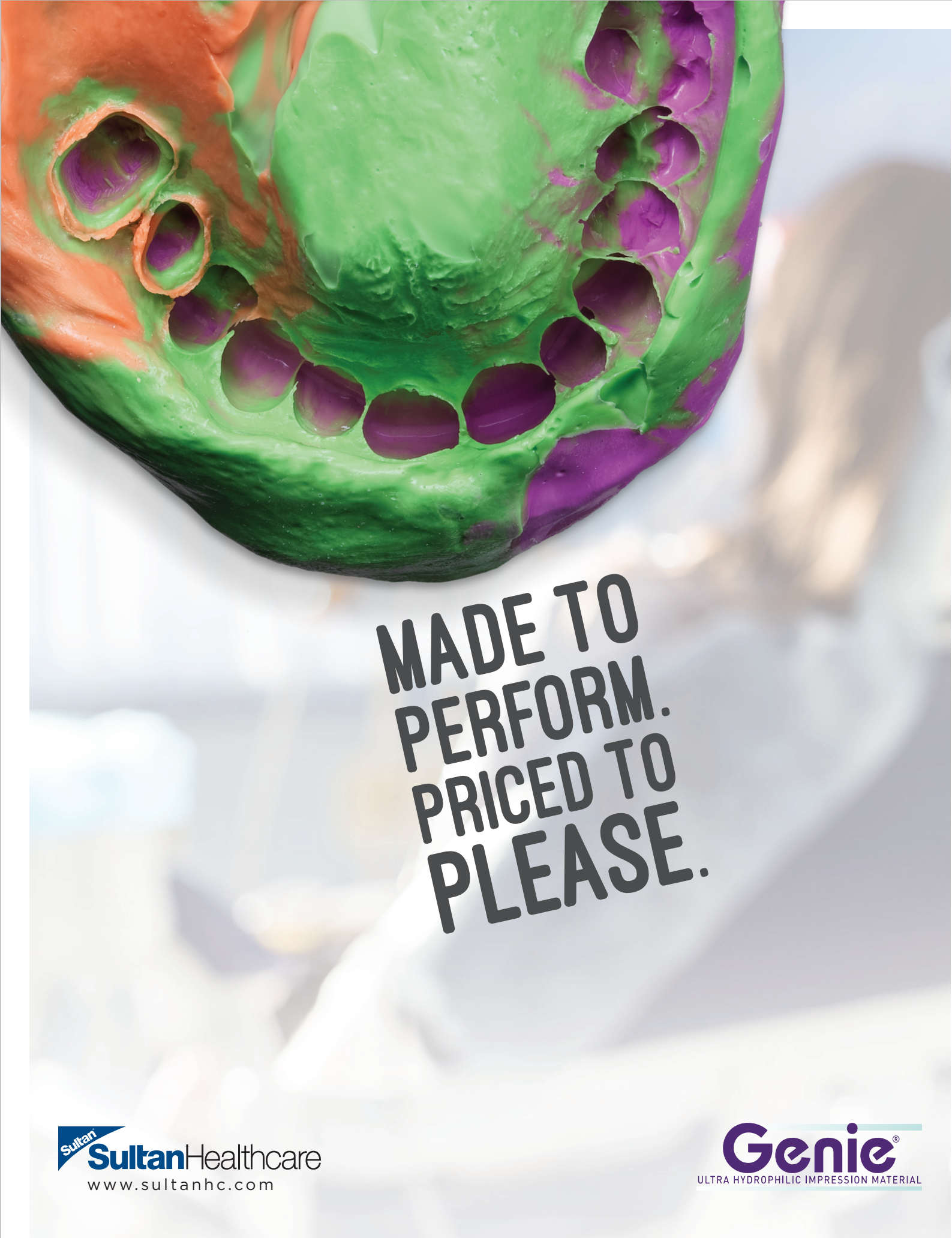

* Based on manufacturers' directions for use.
**Based on compressive strength of leading competitors shown.
†Data on file.
Competitive Brands are not registered trademarks of Sultan Healthcare.



GENIE® MATERIAL OFFERS MORE THAN A GREAT PRICE.

It takes more than a low price to become a top-selling impression material. Genie® material also offers an exceptional experience for dentists and patients, remarkable flowability and reliable performance. If you want more from your impression material, reach for the leading value impression brand.

“Our office uses Genie® material for all of our Orthodontic cases. It tastes good, sets fast and is very accurate.”
— J. Juan, Dental Assistant
SUMMERVILLE, SC

EXCEPTIONAL EXPERIENCES FOR DENTIST AND PATIENT.

Many impression materials can deliver results once. A truly high quality material delivers every time. For 15 years, Genie® material has proven its worth to clinicians with both its performance and its high customer satisfaction. Because we see your practice’s success as the benchmark for our own, Genie material offers a great price and a whole lot more. Want proof? Just listen to what Genie® material users have to say.

- ✓ Choice of five viscosities and two set times (2:10 or 4:30), for excellent flowability and capture of margins
- ✓ Color makes it easy to read
- ✓ Pleasant berry scent
- ✓ Fast set time for patient comfort

93.7%
OF DENTISTS WERE SATISFIED
WITH GENIE MATERIAL’S FLOW⁺

93.0%
OF DENTISTS WERE SATISFIED
WITH GENIE MATERIAL’S
RELIABILITY + PREDICTABILITY⁺

90.2%
OF DENTISTS WERE SATISFIED
WITH GENIE MATERIAL’S
EASE OF HANDLING⁺

Performance Attribute	Overall Dentist Satisfaction ⁺
Price	96.7%
Ability to work in wet environment	94.1%
Material flow	93.7%
Reliability & predictability	93.0%
Material tear strength	92.3%
Broad case suitability	91.8%
Material color/contrast	90.9%
Ability to capture margin detail	90.9%
Ease of handling	90.2%



REMARKABLE FLOWABILITY FOR AMAZING DETAIL.

An impression material should flow past the prepared margin to capture all the detail you and your lab need to deliver quality restorations. Genie® material delivers superior performance over comparably priced competition with up to 2.5 times the flow. Thanks to its excellent flow characteristics under pressure, you can count on the detail you need. Genie® material's impressive performance in the material flow test, shown below, speaks for itself.

“I switched to Genie® material to get better results with my crown and bridge impressions. Genie material provided that improvement in a consistent manner and the lower price was an excellent bonus.”

— S.A. Terry, DDS
North Vernon, IN

Genie® material
26mm

MATERIAL FLOW TEST
A material flow test determines how far material can flow into a narrow mold. Wash material is placed at the base of the mold, then pressure is applied to engage its thixotropic nature. The better the flow, the more complete the resulting fin will be. Testing was conducted at RT & 30 seconds after start of mix. Average of three samples.

3M ESPE Paradigm™
24mm

Patterson Reflection™
20mm

GC Examix™
14mm

DenMat® Splash!® Max
14mm

Clinician's Choice Affinity™
11mm

Kettenbach Panasil®
11mm

GC Exafast™
10mm

Henry Schein VP Mix HP
10mm

DenMat® Splash!®
9mm

Danville Star VPS™
8mm

RELIABLE PERFORMANCE FOR EVERY IMPRESSION.

Just because an impression material is affordable doesn't mean you can't expect it to offer impressive performance in every case. Genie® material offers up to 1.9 times the tear strength and similar hydrophilicity of comparably priced competitors, so you can always count on pulling accurate, dependable impressions. By performing time and time again, Genie® material sets you up for predictability.

- ✓ Excellent flow characteristics under pressure provide better detail and surface coverage
- ✓ Ultra-hydrophilic chemistry designed to work in wet environments
- ✓ Predictable, reliable results case after case
- ✓ Can be used in a variety of cases: standard to multiple implant impressions, single-unit crown, dentures and more.

Sultan Genie® material	Tear Strength TM-000029 (psi)	Contact Angle Uncured Film (2" @ 80% RH)
GC Exafast™	121	68°
3M ESPE Paradigm™	141	73°
Henry Schein VP Mix™ HP (Berry)	132	88°
Clinician's Choice Affinity™	224	78°
Danville Star VPS™	113	68°
DenMat® Splash!®	205	61°

“Genie® material provides reproducibility at a great price, which adds up to value that you can count on.”

— D. Chacko, DMD
ONEIDA, TN



FOLLOW A GREAT GENIE® MATERIAL IMPRESSION WITH A RELIABLE VERSA-TEMP®2 MATERIAL PROVISIONAL.



The temporary is an essential part of the restorative process... important to your patient, and a reflection of the work you and your practice produce. Versa-Temp® 2 material combines the superb performance and natural-looking aesthetics you want—for temporaries that set quickly, handle easily, look great and stay strong until the final crown or bridge arrives.



UP TO
42%
FASTER* SET TIMES THAN
LEADING COMPETITORS

That means less chair time, less expense, and happier patients.

UP TO
70%
STRONGER** THAN
LEADING COMPETITORS

At 450 MPa compressive strength, Versa-Temp® 2 material is formulated to last— to help ensure your patients don't have to come back before their crown or bridge is ready.


EASY TO
SHAPE
WITH BEAUTIFUL RESULTS

Versa-Temp® 2 material makes the temporary process smooth and easy, for aesthetically pleasing temporaries in minutes!

- ✓ Simple to adjust and trim
- ✓ Natural fluorescence
- ✓ Easy to polish with an alcohol wipe (e.g. Isopropyl)



**4 MINUTE
TOTAL SET TIME!**

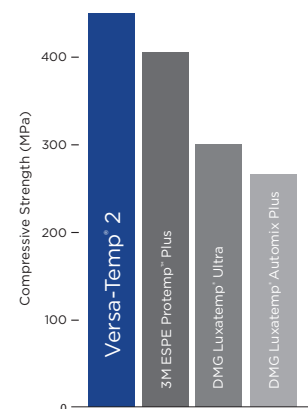

A GREAT PROVISIONAL BEGINS WITH SILGIMIX® MATERIAL.



An alginate replacement combines the low cost of an alginate with the detail and versatility of a VPS. Silgimix® material works hand-in-hand with Versa-Temp® 2 material. The result is a great impression with less time adjusting and polishing... and a high-quality temporary.

THE ECONOMICAL AND RELIABLE SULTAN HEALTHCARE FAMILY OF PRODUCTS.



RESTORATIVE MATERIALS



TEMPORARY CROWN AND
BRIDGE MATERIALS